

# Speedex putty soft

## Instructions for use

EN

Please read the Instructions for Use carefully before using the product. Keep for later reference.

### PRODUCT DESCRIPTION

Speedex is a two-component impression material based on condensation-type polysiloxanes. After manual mixing of base and activator, the Speedex materials form pastes which are used individually or in combination as dental impression material, usually together with a standard commercial or individual impression tray and using conventional impression techniques.

### INTENDED USE

Speedex is intended for the recording of the current physical situation in a patient's mouth („dental impression material“) for the purpose of repairing, reshaping or replacing the patient's teeth.

### COMPOSITION

Condensation-type silicone elastomer, putty material.

- Polysiloxane
- Colour: light blue

### INDICATION

- Primary impression in the putty-wash impression technique
- Pick-up impression in the simultaneous mixing technique
- Impressions for study models, orthodontic models, keys and protective isolation of teeth during fitting of prostheses

### CONTRAINDICATION

- Do not use in case of known allergies to ingredients of the Speedex impression materials.
- Loose teeth can be further loosened or extracted by taking an impression
- Material combinations with impression materials of other manufacturers are not permitted

### SAFETY INSTRUCTIONS

- Eyes, respiratory organs and skin can be irritated by Speedex Universal Activator
- Impressions that came into contact with the oral mucosa may be contaminated. Observe disinfection instructions
- Only supplied to dentists and dental laboratories or upon their instructions
- Contains nano materials (bonded particles)

### SIDE EFFECTS / INTERACTIONS

No harmful reactions or secondary effects on patients and / or dental personnel are known.

### USER- / PATIENT GROUP

The product shall be used by qualified dental professionals only. Suitable for all patient groups including children, elderly people and pregnant women.

### PREPARATION


#### Trays

All impression trays that have been designed for dental use, either prefabricated or individual, can be used for impression taking. For perfect adhesion, we recommend applying a thin layer of COLTENE Adhesive to the trays before use. This prevents detachment of the impression material during removal from the oral cavity. Detailed information on the application can be found in the COLTENE Adhesive Instructions for Use.

### PROPER USE

#### Dosage

Dose Speedex putty soft with enclosed measuring-scoop (scrape off surplus). Spread out the material in the palm of hand. For each measuring scoop press the scoop rim once into the material on the hand. Apply one strand length of Speedex Universal Activator per circle diameter. Make sure that the tins and tubes are closed immediately after use.

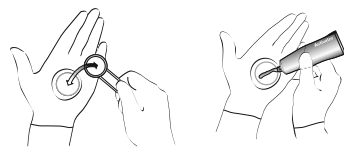

⚠ Take the following safety precautions when working with Speedex Universal Activator:

- Protective glasses are recommended. In case of contact with eyes, rinse immediately with water and consult an ophthalmologist
- Avoid direct skin contact
- Always wear protective gloves when working (vinyl gloves are recommended)
- Avoid inhalation or ingestion

### APPLICATION TIME

#### Normal dosage



Intensive kneading, high temperatures and humidity accelerate the curing of the impression materials. Low temperatures slow down the curing process. By over- and underdosing the Speedex Universal Activator the setting time can be accelerated or slowed down (max. ± 20% of the normal dosage)

#### 20% Activator overdosage



#### 20% Activator underdosage



### Mixing

Fold material and knead with fingers until an absolutely homogeneous colour is attained, but for at least 30-45 seconds mixing time. The use of a stopwatch is recommended for this purpose. When using mechanical mixing devices, observe the instructions of the device manufacturer. After successful (streak-free) mixing, form the material into a strand and place it immediately into the impression tray.

### Impression

Insert the impression tray, loaded with impression material, into the oral cavity and press for 2-3 seconds. Then hold in situ until the material has set completely.

Avoid using excessive amount of silicone for impression taking due to the risk of overflow material being swallowed.

### TESTMETHOD FOR CORRECT APPLICATION

Material residue needs considerably more time to set at room temperature than in the patient's mouth. You should therefore check intraorally if the material is set before removing from the mouth. If the manual pressure test does not leave a pressure mark in the impression material after the time specified in the chapter on „Dosage“, the impression can be removed from the patient's mouth. However, if the test leaves a visible pressure mark in the material, the material is not yet fully cured.

### REPROCESSING, CLEANING, DISINFECTION AND MAINTENANCE

#### Cleaning of the Impression

After the putty impression has been removed from the patient's mouth, it must be thoroughly cleaned before further use with wash materials (Speedex light body/medium). After cleaning, the final impressions can be disinfected with commercially available disinfectant solutions (e.g. Dürr Dental SE, MD 520 impression disinfection) specifically designed for dental impression materials and dried using air pressure. If H<sub>2</sub>O<sub>2</sub> hydrogen peroxide is used for disinfection, it is recommended to thoroughly rinse with lukewarm water afterwards in order to avoid bubble formation.

#### Cleaning of trays

The impression can be removed with a blunt instrument. Adhesive residues can be dissolved by placing in commercially available universal solvents or benzine.

#### Manufacture of models

The ideal time lies between 30 min and 72 h after taking the impression. The surface tension will be reduced and pouring will be facilitated if the impression is briefly washed out with a detergent and rinsed thoroughly in lukewarm clear water afterwards. All industry-standard dental stone model materials (i.e. Fujirock Dental Stone, Hard Rock Dental Stone) can be used.

#### SHELF LIFE/STORAGE

- Expiry date: see label on primary packaging
- Storage temperature: 15-23°C / 59-73°F
- Relative humidity: 50% ± 10%
- Shelf life after first opening, 3 months

Protect against exposure to heat and sun. Avoid extreme temperature fluctuations. Close tins and tubes immediately after use.

### DISPOSAL

Dispose of waste according to applicable legislation. Special country-specific regulations may apply. Can be disposed together with household waste in compliance with official regulations in contact with approved waste disposal companies and with authorities in charge. (Only dispose of completely emptied packages.)

### TECHNICAL DATA

#### Technical data according to ISO 4823

The technical data were determined at 23°C/73°F and 50% relative humidity.

Mixing time:	00:30 min
Processing time:	≤1:15 min
Oral setting time:	3:00 min
Mixing ratio Putty / Activator:	23 ml / 1 ml
	resp. 6,90 g / 0,15 g

### REPORTING OBLIGATION

All serious incidents occurring in conjunction with this product must be reported immediately to the manufacturer as well as to the competent authority.

In the unlikely event of inhalation, ingestion, eye contact, or similar incidents seek immediate medical attention from an appropriate medical specialist to mitigate potential harm.

### SAFETY DATA SHEET

www.coltene.com



GLOSSARY	
	Consult instructions for use
	Keep away from sun light
	Keep dry
	Temperature limitation
	Marking of Conformity Europe
	Identification for Ukraine
	Legal Manufacturer
	Expiry Date
	Batch Code
	RX only
	Medical Device
	Manufacturing Date
	Unique Device Identifier
	European Authorized Representative
	Importer
	Reference Number

### COLTENE International Dental Group

Dent4You AG  
Bahnhofstrasse 2  
CH-9435 Heerbrugg

Manufactured by  
Coltene/Whaledent AG  
Feldwiesenstrasse 20  
CH-9450 Altstätten

Customer Center  
service@coltene.com

**COLTENE**