

3 GALLON / 11.4 L ULTRASONIC CLEANING SYSTEM

BioSonic[®] UC300 / UC300R

Your Patients Trust You
You Can Trust **BioSonic**



 **COLTENE**

BioSonic® UC300/300R




FLEXIBLE ORGANIZATION

- Multiple basket and beaker arrangements
- Adjustable basket dividers for instrument sorting
- Accessory basket accommodates popular cassette systems



MAXIMUM PERFORMANCE

- Cover traps aerosols and directs condensation back into tank
- Externally accessible drain allows for easy emptying
- Quiet operation for a comfortable work environment

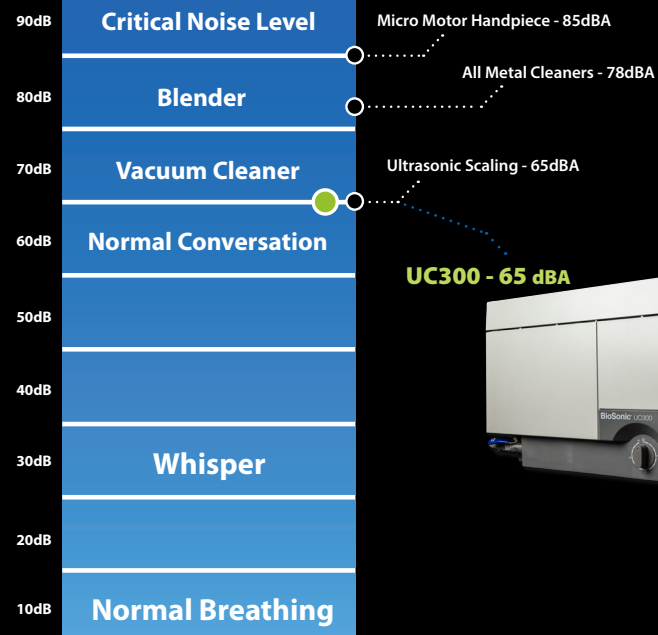

TRUSTED TECHNOLOGY

Cavitation is the heart of ultrasonic cleaning technology. This occurs as millions of bubbles, generated from high frequency soundwaves, form and burst continuously in the liquid. When the bubbles collapse, energy is released in the form of heat and pressure waves, thus breaking down the bonds that hold heavy debris on your instruments. The combination of tiny pressure waves created by cavitation and enhanced by specially formulated solutions reach deep into crevices, removing gross and microscopic debris. The BioSonic UC300 ensures safe and reliable device operation for your daily infection control needs.







VERSATILE CONFIGURATION

- Generous 3 gallon / 11.4 liter use capacity
- Countertop and recessed models available
- Optional remote digital timer for user-friendly operation on recessed model



BioSonic® UC300/UC300R Ultrasonic Cleaning System

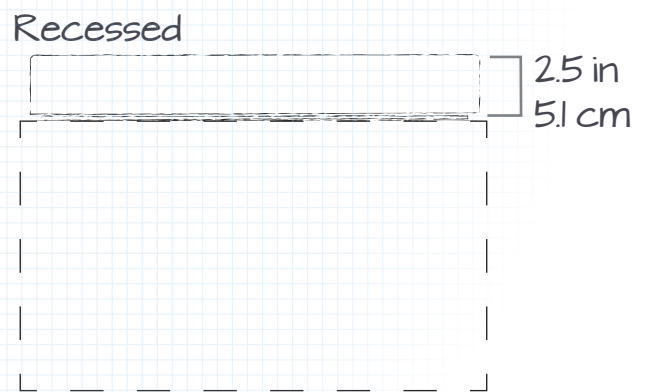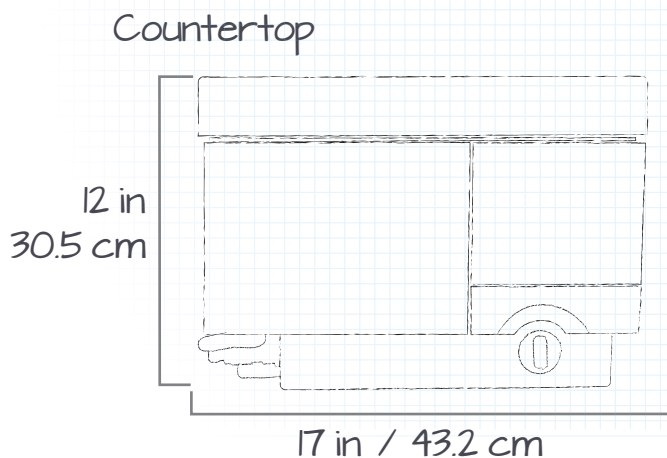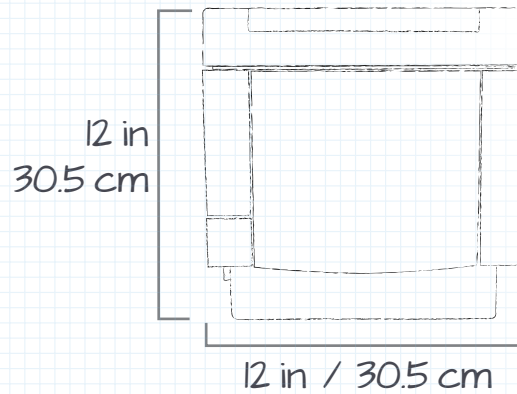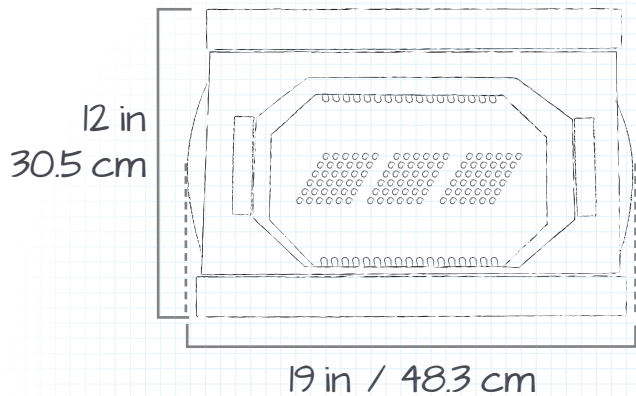

For maximum performance and longevity of your BioSonic Ultrasonic Cleaner, we recommend the use of genuine BioSonic solutions and accessories.

Order Information

UC300B115	BioSonic Ultrasonic Cleaner Accessories not included
UC300RB115	BioSonic Recessed Ultrasonic Cleaner Accessories not included
UC310	Large FingerGuard™ Instrument Basket
UC311	Twin FingerGuard™ Instrument Basket
UCB2	Small Parts Basket
UC315115	Remote Digital Timer for UC300R
UC325	Cassette Basket for UC300 / UC300R
UC330	Two-Beaker Holder
UC53	Borosilicate Beaker 600 ml with cover and positioning ring
UC335	Basket Section Inserts

BioSonic Solutions

UC30 General Purpose Cleaning Solution - 16 oz	Contaminants: Blood, debris, saliva and general cleaning	Items: Instruments, crowns, bridges, dentures and ultrasonic scaler tips
UC32 Enzymatic Ultrasonic Cleaning Concentrate - 8 oz	Contaminants: Protein soils (dried blood, tissue, saliva) and general cleaning. Effectively breaks down and loosens hardened debris	Items: Instruments, crowns, bridges, dentures and ultrasonic scaler tips
UC34 Plaster and Stone Remover Powder - 1.85 lb	Contaminants: Plaster and stone	Items: Dentures, partials, castings, denture flasks, mixing bowls and non-aluminum articulators
UC39 Ultrasonic Cement Remover Concentrate - 16 oz	Contaminants: Temporary cements and permanent cements (including zinc-phosphate, polycarboxylate)	Items: Crowns, veneers, brackets and instruments
UC42 Germicidal Ultrasonic Cleaning Concentrate - 16 oz	Contaminants: Blood, debris, saliva and general cleaning. Cleans instruments as it disinfects the ultrasonic bath	Items: Instruments, crowns, bridges (no plastic parts)