

BRILLIANT EverGlow® Flow



THE IDEAL COMPLEMENT

BRILLIANT EverGlow Flow is based on the same filler technology as BRILLIANT EverGlow and is therefore the perfect complement. Thanks to its well-adjusted thixotropy, BRILLIANT EverGlow Flow has a low viscosity yet offers high stability. A new, extra thin application needle extends the range of accessories and allows for the highest application precision.

FUNCTIONAL SHADE CONCEPT – INCLUDING OPAQUE SHADE

BRILLIANT EverGlow Flow comes with all the quality characteristics of BRILLIANT EverGlow filler technology, such as polishability and gloss retention. The shade assortment has also been adopted and offers the convenient Duo Shades and three translucency levels.




Although the universal and translucent shades serve most cases, dark areas and chroma deviations can negatively affect the final aesthetic outcome of a restoration and demand for an opaque cavity lining. The Opaque A2 shade was developed to suit this need by creating the basis for a harmoniously shaded and aesthetic restoration.

As a result, BRILLIANT EverGlow Flow delivers convincing aesthetic outcomes and is the perfect partner for the universal composite BRILLIANT EverGlow.

HIGHLY STABLE, YET VERY FLOWABLE

With BRILLIANT EverGlow Flow, excellent stability and good flowability are no longer a contradiction. BRILLIANT EverGlow Flow manages to combine these properties in a single material which makes it extremely versatile to use. With its high stability and optimal thixotropic property, BRILLIANT EverGlow Flow can be applied with absolute precision in a controlled and easy manner.

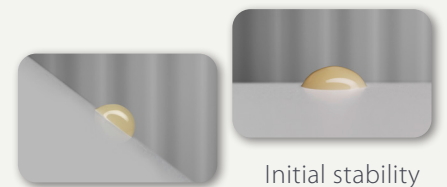

GOOD WETTABILITY FOR EXTENSIVE COVERAGE

While the controlled flow behavior facilitates focused work, the low viscosity of the material provides good wetting which, for example, is important in the case of cavity linings and small fillings.



TWO APPLICATION NEEDLE SIZES FOR MULTIPLE INDICATIONS

Every BRILLIANT EverGlow Flow syringe comes with two application needle sizes. In addition to the standard size with an inner diameter of 0.8 mm, the accessories provided have been extended with a fine size. With only 0.4 mm inner diameter, it allows for the highest application precision.



ORDER INFORMATION

BRILLIANT EverGlow Flow Intro Kit	Opaque A2 (OA2)	Bleach (BL)	A1 / B1	A2 / B2	A3 / D3	A3.5 / B3	A4 / C4	Translucent (Trans)
60019760 4 x 2g syringes			1	1	1	1		

Kit accessories: Instructions for use, mixing pad, 40 application needles

REFILLS

BRILLIANT EverGlow Flow Refills (1 x 2 g + 10 application needles)

60019751	Opaque A2 (OA2)	60019755	A3 / D3
60019752	Bleach (BL)	60019756	A3.5 / B3
60019753	A1 / B1	60019757	A4 / C4
60019754	A2 / B2	60019758	Translucent (Trans)

8469	20 Application needles Standard (Ø 0.8 mm)
60022137	20 Application needles Fine (Ø 0.4 mm)



COMPLEMENTARY PRODUCTS

ONE COAT 7 UNIVERSAL

Light curing single-component bond for adhesive luting with the self-etch, total etch or selective etch technique.

60019539 ONE COAT 7 UNIVERSAL Refill, 5 ml



Alpen ShapeGuard® Composite Plus

Silicone polisher for pre-polishing (purple) and high gloss polishing (blue) of composite restorations. Due to its flexible lamellae, the ALPEN ShapeGuard adapts to the surface structure.

60021910 Alpen ShapeGuard® Composite Plus Assorted Pack



Alpen ShapeGuard® Composite Plus

Coltène/Whaledent Inc.
235 Ascot Parkway
Cuyahoga Falls, Ohio 44223 / USA
T +1 330 916 8800
F +1 330 645 8704
info.us@coltene.com

Coltène/Whaledent AG
Feldwiesenstr. 20
9450 Altstätten / Switzerland
T +41 71 757 53 00
F +41 71 757 53 01
info.ch@coltene.com

www.coltene.com © 2020

COLTENE

P/N 31626A