


ONE-COMPONENT LIGHT-CURED UNIVERSAL ADHESIVE

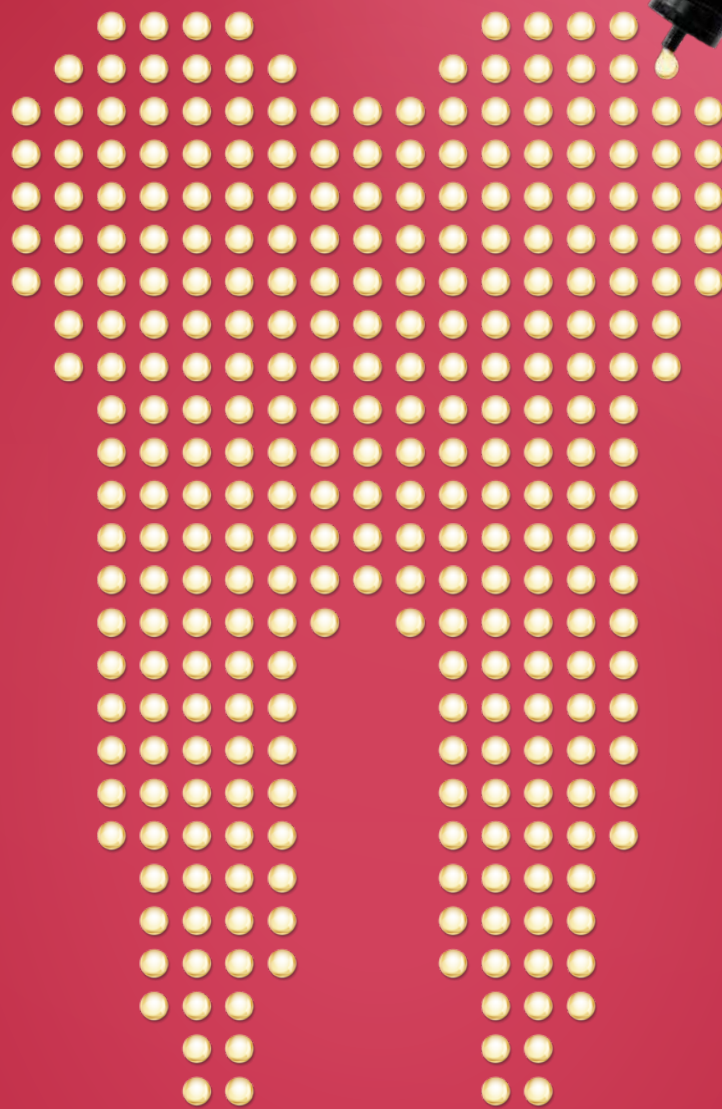

# ONE COAT 7 UNIVERSAL

The bond for all situations

 COLTENE

The light-curing **ONE COAT 7 UNIVERSAL** is the only bond you need for direct and indirect filling therapy in your practice. No priming, no silanising – just apply, disperse with air and cure – finished.

Advantages at a glance



Only one universal bond for all indications and etching techniques



Extremely economical: up to 300 applications per bottle




Exceptional bonding values and excellent marginal seal

The all-rounder


# ONE COAT 7 UNIVERSAL

Reliable bonding on different surfaces

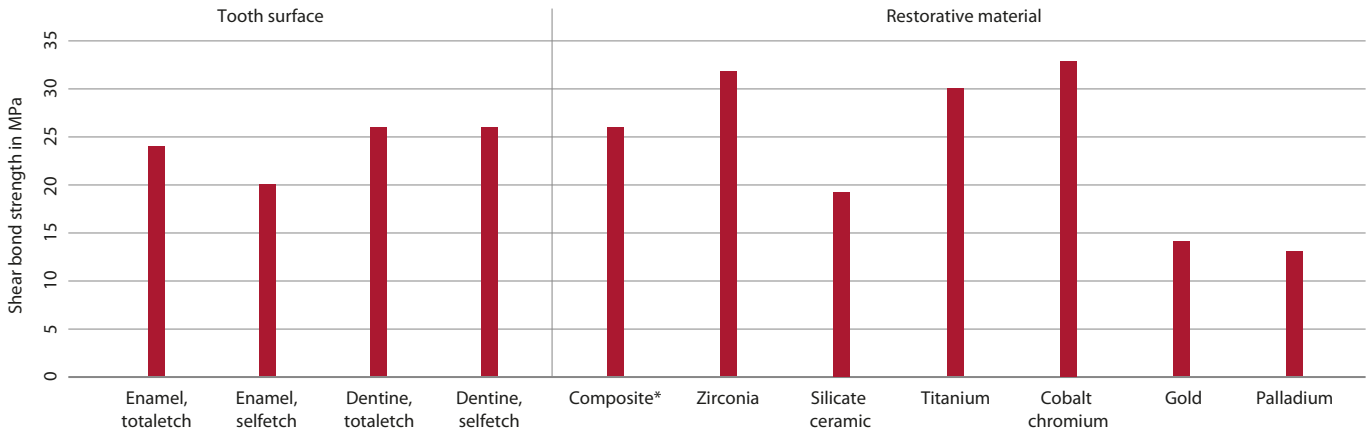

\*SYNERGY D6, COLTENE  
Source: internal data  
Measurement according to Watanabe

Tooth substance, ceramic, metal, composite – ONE COAT 7 UNIVERSAL gives you safe and reliable adhesion on any surface. Depending on the situation, you can either choose to use the self-etching effect or additionally etch enamel and / or dentine with phosphoric acid (e.g. Etchant Gel S, COLTENE). Are you luting an indirect restoration and require a chemically curing composite material for this purpose? No problem – the available Activator ensures compatibility with ONE COAT 7 UNIVERSAL.

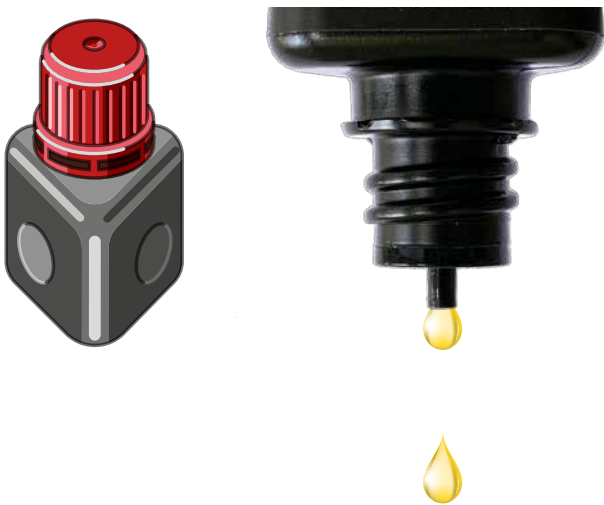

Optimised application for comfortable working and economical use

The ergonomic triangular bottle sits perfectly in the hand. The optimised dropper impresses with precise, clean dosing as well as application without dripping. This allows up to 300 applications per bottle.

The perfect bonding partner  
for the COLTENE range.



## Order Information

### ONE COAT 7 UNIVERSAL

#### 60019538 Intro Kit

1 × 5 ml	ONE COAT 7 UNIVERSAL
1 × 2.5 ml	Etchant Gel S
1 × 20	Application tips
1 × 50	Brushes black
1 ×	Mixing well
1 ×	Cover for mixing well

#### 60019539 Refill

1 × 5 ml	Refill
----------	--------

#### 60019540 Single Dose

50 × 0.1 ml	Single Dose
1 × 50	Brushes black

#### 7634 Brushes black, 50 pcs.

#### 7054

1 × 3 ml	One Coat 7.0 Activator
1 × 50	Brushes black

#### 7373

2 × 2.5 ml	Etchant Gel S Kit
1 × 20	Syringe
	Application tips

© COLTENE – [www.coltene.com](http://www.coltene.com)



#### Coltène/Whaledent Ltd.

The President Suite  
Kendal House, Victoria Way  
Burgess Hill, West Sussex  
RH15 9NF / UK  
T +44 1444 235486  
[info.uk@coltene.com](mailto:info.uk@coltene.com)

#### Coltène/Whaledent AG

Feldwiesenstrasse 20  
9450 Altstätten / Switzerland  
T +41 71 757 53 00  
[info.ch@coltene.com](mailto:info.ch@coltene.com)

 **COLTENE**