



**DuoTEMP<sup>®</sup>**   
*dual cure temporary filling material*

## Double security keeps you smiling!

- Dual-cure
- High marginal integrity
- Easy removal
- One material with multiple indications

**DuoTEMP is a dual-cured, eugenol-free, white temporary filling material formulated from zinc oxide/zinc sulphate for optimal protection from the infiltration of bacteria. DuoTEMP is one of the first all-in-one temporary filling materials that can be cured using both light and saliva. This solution delivers greater flexibility, allowing the patient to immediately resume normal daily activities. DuoTEMP also contains a pleasant flavor that keeps the dentist and patient smiling.**

► **Double security**

• **Dual-cure**

After only 40 s of light curing, the patient can immediately load the material. Afterwards, DuoTemp will continue to completely self-cure on its own with saliva.

• **High marginal seal due to slight expansion**

Reduces infiltration of bacteria

► **One material with multiple indications**

Reduces inventory and minimizes confusion

- Temporary fillings of class I and II cavities
- Temporary fillings for endodontic applications
- Temporary fillings for inlay / onlay preparations

► **Natural tooth-like feeling**

DuoTEMP's high surface hardness feels like a natural tooth to minimize irritation and maximize patient comfort.

► **Additional benefits**

- Contains zinc oxide and fluoride for additional protection
- Eugenol-free – no effect on the curing of composites
- Radiopaque

**Order Information:**

**5830** DuoTEMP Intro Kit

2 syringes of 5 g, Intraoral Separator 5ml

**5831** DuoTEMP Single Pack


1 syringe of 5 g

**7090** DuoTEMP Eco Pack

5 syringes of 5 g

**7235** Intraoral Separator, 5ml

► **Marginal Seal Comparison**

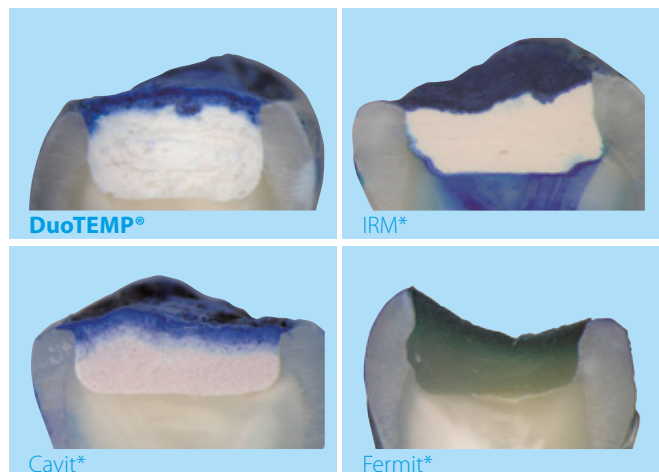

Cross-sectioned extracted human molars were filled with cured DuoTEMP and then immersed in methylene blue dye to determine the marginal sealing properties. DuoTEMP demonstrated excellent sealing properties.

Source: internal data

\* Not a trademark of Coltène/Whaledent

► **Easy, time saving application**



One component without mixing. Good retention in the cavity. Easy removal.

© 2015 Coltène/Whaledent AG – www.coltene.com

**Coltène/Whaledent AG**

Feldwiesenstrasse 20  
9450 Altstätten / Switzerland  
Tel +41 71 757 5300  
Fax +41 71 757 5301  
info.ch@coltene.com

**Coltène/Whaledent Inc.**

235 Ascot Parkway  
Cuyahoga Falls, Ohio 44223 / USA  
Tel +1 330 916 8800  
Fax +1 330 916 7077  
info.us@coltene.com

**Coltène/Whaledent Ltd.**

The President Suite-A  
Kendal House, Burgess Hill  
Victoria Way  
West Sussex, RH15 9NF/U.K.  
Tel. +44 (0)1444 235486  
Fax +44 (0)1444 870640  
info.uk@coltene.com

