



CanalPro[®]
Endodontic Irrigation System



COLTENE Covers the Entire Endo Procedure

COLTENE offers a complete line of endodontic products for the root canal procedure, from diagnosis through restoration.

Visit coltene.com/products/endodontics to view the entire endo line.

CanalPro®

Endodontic Irrigating Solutions



Saves Times and Improves Outcomes

When it comes to irrigating, CanalPro Endodontic Solutions are your best choice for successful treatment. We have engineered our endodontic solutions to minimize the time spent on irrigation, giving you the best approach to cleansing canals and the best outcomes.



Innovative Delivery System

Innovative, modular system for the safe and efficient filling of syringes with irrigation solutions.



Color-coded for Safety and Convenience

Organizing procedures and products by color generates a more efficient work environment.

CanalPro NaOCL

- For irrigation/debridement of root canals during and after instrumentation
- Available in 6% and 3% formula



CanalPro EDTA

- 17% EDTA solution (pH 8.5)
- Removes smear layer and dentin mud
- Opens dentin tubules for:
 - Disinfection solutions to work better
 - Better adhesion of sealers and obturation materials



CanalPro CHX-Ultra

- Same patent pending formula as CHX-Plus™*
- Use as a final step for long-lasting cleansing
- Kills multi-species biofilm two times faster and planktonic bacteria ten times faster than standard 2% Chlorhexidine
- Powerful wetting agents and surface modifiers enable CanalPro CHX-Ultra to penetrate into lateral canals and isthmuses
 - Multiple Patents Pending



Catalog Number	Description	Package Size
60011161	CanalPro NaOCL	6%, 480 ml
60011160	CanalPro NaOCL	3%, 480 ml
60011157	CanalPro EDTA	480 ml
60011158	CanalPro EDTA	120 ml
60013994	CanalPro CHX-Ultra	480 ml

Color Syringes

- Helps increase safety and minimize the chance of syringe swap
- Latex-Free, color-coded syringes offer a fast, easy way to organize and identify syringes for irrigants and solutions
- Standard luer lock style



Catalog Number	Description	Package Size
60011173	Red, 10 ml	50 syringes
60011174	Blue, 10 ml	50 syringes
60011175	Yellow, 10 ml	50 syringes
60011176	White, 10 ml	50 syringes

Catalog Number	Description	Package Size
60019321	Red, 5 ml	50 syringes
60019322	Blue, 5 ml	50 syringes
60019323	Yellow, 5 ml	50 syringes
60019324	White, 5 ml	50 syringes

CanalPro Syringe Station

- No Mess! No Waste! Less Time!
- Cost-effect way to create a dedicated area for syringe filling
- Syringe activated valve
- Low level light indicator
- Helps prevent counter, carpet and scrub damage
- Several units may be connected to create a convenient station for filling multiple solutions



Catalog Number	Description
60011178	CanalPro Syringe Station

CanalPro Flex Tips

Bendable to facilitate easy access for irrigation, application and micro-aspiration.


- Excellent canal tracking
- Unique non-kink polyamide tip
- Disposable (single use) tip
 - Patented



Catalog Number	Description
	CanalPro Flex Tips
60011164	30 ga 20 per pack

CanalPro Slotted-End Irrigation Tips

- Ideal for irrigation of canals, pockets and fistulas
- Slotted and side-vented

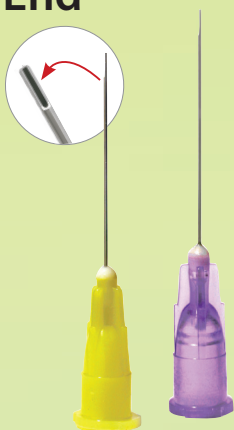

Catalog Number	Description	
	CanalPro Slotted-end Irrigation Tips	
60011167	27 ga	100 per pack
60011168	30 ga	100 per pack

CanalPro Side-Port Irrigation Tips

- Ideal for irrigation of canals, pockets and fistulas
- Closed end for side port delivery




Catalog Number	Description	
	CanalPro Side-port Irrigation Tips	
60011169	27 ga	100 per pack
60011170	30 ga	100 per pack

SurgiTip-endo®

Endodontic aspirator tip used to quickly dry root canals.

- Unique, multi-directional tip design offers flexibility without kinking
- Simply place into saliva ejector, insert into canal and aspirate
- After use, only 1-2 paper points are needed to dry canal
- Supplied sterile, ready to use

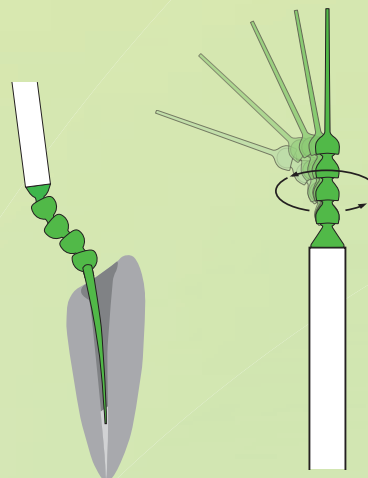

Catalog Number	Description	Package Size
60018280	0.35 mm tip opening 6.5 mm end opening	20